





sartorius stedim
biotech

Vivapure® Ion Exchange Chromatography Products

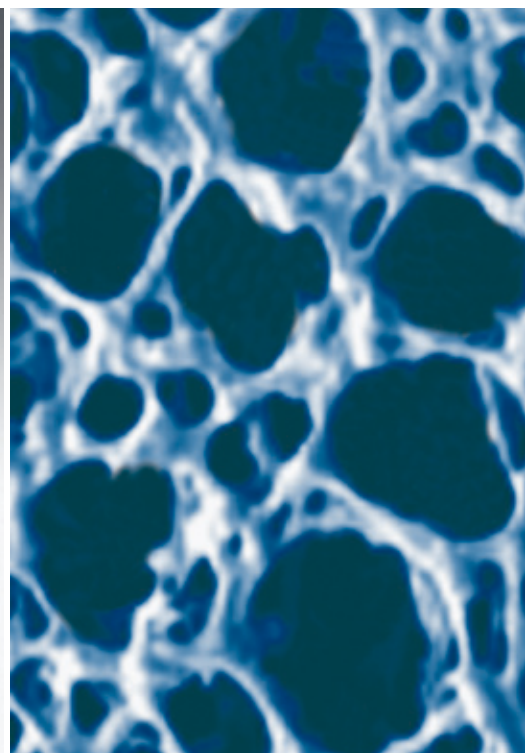

turning science into solutions

Vivapure® Ion Exchange Protein Purification Products

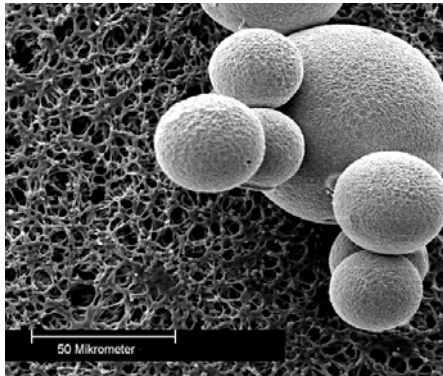

Chromatography gel beads (right) are shown on top of a membrane adsorber in this SEM picture. The membrane adsorber pores are over 50 × larger than bead pores.

The unique Sartobind membrane adsorber matrix

Sartobind IEX membrane adsorbers are based on stabilized regenerated cellulose and display a microporous structure with a pore size of > 3 µm, which is orders of magnitude larger than conventional chromatographic gel materials. This allows molecules to be transported to the ligands immobilized on the membrane adsorber by convective flow, leading to very high flow rates.

In contrast to that, gel chromatography is slowed down due to diffusion limitations, as the molecules need to enter the small bead pores in order to be bound by the ligands. The porous membrane adsorber enables fast, reproducible and scalable protein purification.

Fast and easy-to-use chromatography devices

Vivapure® Ion Exchange (IEX) products, incorporate Sartobind Membrane Adsorber technology as their chromatography matrix. They make protein purification as easy as filtration. The devices are ready-to-use and do not bear the risk of running dry. For many protein purification applications, they can replace time-consuming and tedious column chromatography.

The rapid 1-2-3 bind-wash-elute protocol especially lends itself to screening applications, where many different samples are processed in parallel.

... Scalable to Industrial Applications



Scalable to industrial scale applications

Sartobind membrane adsorber (MA) technology is available from 96-well plate format to industrial process scale modules. Sartobind IEX MA units are useful tools when flow control is needed, as the down scale step to a process scale application.

Scalable SingleSep disposable capsules for process scale applications, e.g. virus removal, DNA removal and endotoxin removal. The devices are available in formats ranging from 1 ml to 1620 ml membrane volume, for capacities from 0.029 g to 48 g. For technical data and description, refer for SingleSep brochures and data sheets.

Vivapure® IEX products



Vivapure® Mini – 400|500 µl
Binding capacities: 1–4 mg



Vivapure® Maxi – 19|20 ml
Binding capacities: 15–80 mg

Membrane availability

Functional groups	Ion exchanger type
Sulphonic acid (S)	Strong acidic cation exchanger: $\text{R-CH}_2\text{-SO}_3^-\text{Na}^+$
Quaternary ammonium (Q)	Strong basic anion exchanger: $\text{R-CH}_2\text{-N}^+(\text{CH}_3)_3\text{Cl}^-$
Diethylamine (D)	Weak basic anion exchanger: $\text{R-CH}_2\text{-NH}^+(\text{CH}_2\text{H}_5)_2$

Available formats

Vivapure®	Application
Vivapure Mini Spin Columns	<ul style="list-style-type: none">– Sample fractionation– Purification condition scouting– Small scale purification
Vivapure Maxi Spin Columns	<ul style="list-style-type: none">– Large scale sample fractionation– One step protein purification concentration– Polishing of his-tagged protein

Vivapure® advantages

Fast and simple to use

- Devices are ready to use – no column packing
- Make protein purification as simple as filtration

Reproducible results

- Membrane adsorber columns cannot crack or run dry
- Membrane adsorber columns are highly reproducible to manufacture

Centrifugal devices

- Offer the possibility of working in parallel

Low bed volume

- Small membrane adsorber bed volumes allow working with lower buffer amounts, leading to concentrated elution fractions

Up-scalable product range

- Process scale modules are available with the same Sartobind IEX membrane adsorber matrix

Applications

- Fractionation of protein mixtures prior to 2D-PAGE
- Purification scouting of an unknown protein
- Sample preparation prior to 1-D or 2-D PAGE
- Removal of endotoxins from monoclonal antibodies
- Preparation of heme moiety from heme containing protein for functional analysis
- Polishing of His-tagged proteins
- General protein purification
- Detergent removal from protein solutions
- Purification of antibodies from serum, ascites or tissue culture supernatant
- Establishing a purification | prepurification protocol for a given protein
- HPLC|FPLC sample preparation
- Purification of membrane-bound proteins



Technical Specifications

Performance characteristics

Device	Protein binding capacity* (mg)	Max. volume per centrifuge run using a swing-out rotor (ml)	Max. volume per centrifuge using a fixed angle rotor run (ml)
Vivapure Mini H	4	0.4	
Vivapure Maxi H	60–80	19	10.5

* Actual yields depend on specific protein sample and selected pH and salt conditions. Yields established using 1 mg/ml BSA in 25 mM Tris/HCL pH 8.0 with Vivapure Q & D spin columns and 1 mg/ml cytochrome c in 25 mM sodium acetate buffer pH 5.5 with Vivapure® S & C spin columns.

Capacities and dimensions

Device	Bed Volume (µl)	Membrane Area (cm ²)
Vivapure Mini H	240	7.48
Vivapure Maxi H	2700	84.40

Membrane Adsorber

Nominal pore size	3–5 µm (Large pore size prevents gel filtration effects and minimizes non-specific adsorption)
Thickness	230–320 µm
Amount of ionic groups (µEquivalents/ml)	145–218 µEquivalents/ml for monovalent ions (D, Q & S) 72 µEquivalents/ml for monovalent ions (C)
Working pH (D & C)	4–10
Working pH (Q & S)	2–12
Approximate pka of ionic groups	D–9.5 Q–11 S–1 C–4.5

Materials of construction

Device	Polypropylene
Supporting matrix	Stabilized regenerated cellulose

Ordering Information

Cat Number	Vivapure Mini Spin Columns	Spin Centrifuge	Columns Tubes
VS-IX01SQ16	Vivapure Mini S & Q H Starter Kit (8 of each ion exchange class)	16	32
VS-IX01DH24	Vivapure D Mini H	24	48
VS-IX01QH24	Vivapure Q Mini H	24	48
VS-IX01SH24	Vivapure S Mini H	24	48

Cat Number	Vivapure Maxi Spin Columns	Spin Centrifuge	Columns Tubes
VS-IX20QH08	Vivapure Q Maxi H	8	16
VS-IX20SH08	Vivapure S Maxi H	8	16

Vivapure IEX Spin Columns – Applications

1. Fractionation of complex protein lysates with IEX-membrane spin columns improves resolution of 2-D PAGE.

Perkin Elmer, Boston (USA) Mary Lopez et al., American Biotechnology Laboratory

2. A fast and simple protocol for finding out the optimal purification conditions (purification scouting) of an unknown protein, using different buffers and IEX spin column chemistries in parallel.

Vivascience, Hanover, Germany C. Neumann et. Al.

3. A simple, fast and reliable method using Vivapure Q Mini for removing highly charged contaminant from samples prior to 2D-PAGE e.g. proteoglycans from cartilage explant.

Roche, Palo Alto, USA

4. Vivapure centrifugal anion exchangers were used to remove endotoxin from research grade monoclonal antibody solutions easily with high protein recovery.

Cambridge Antibody Technology, Cambridge UK B. Fish et al., Drug Discovery 3 (2003) 26-27

5. A simple and reliable protocol for the separation of Tween 20 treated soluble Guanylate cyclase from detergent and free porphyrin is described using Vivapure Mini Q for routine use.

Bayer Ag Wuppertal, Germany P. Schmidt et al., Protein Expression and Purification 31 (2003) 42-46

6. Application of Vivapure S centrifugal ion-exchange membrane devices for the purification/polishing of Histagged proteins for crystallization.

Birkbeck College, University of London, UK.

7. Purification of caspase-14 using Vivapure Mini H D anion exchange spin columns for functional analysis.

Department of Medicine, University of Washington, Seattle, USA Andy J. Chien et al., Biochem. And Biophys. Research Communications 296 (2002) 911-9177.

8. Rapid removal of the detergent, n-octyl beta-D-glucopyranoside from a membrane protein mimic using an innovative centrifugal anion exchange membrane technology.

University of Cambridge, UK.

9. Purification of monoclonal antibodies from high cell density (MiniPERM) tissue culture supernatant using Vivapure centrifugal ion exchange membrane devices.

University of Bochum, Germany.

Available in detail on our website:
www.sartorius-stedim.com



Sales and Service Contacts

For further contacts, visit www.sartorius-stedim.com

Europe

Germany

Sartorius Stedim Biotech GmbH
August-Spindler-Strasse 11
37079 Goettingen

Phone +49.551.308.0
Fax +49.551.308.3289

www.sartorius-stedim.com

Sartorius Stedim Systems GmbH
Schwarzenberger Weg 73-79
34212 Melsungen

Phone +49.5661.71.3400
Fax +49.5661.71.3702

www.sartorius-stedim.com

France

Sartorius Stedim Biotech S.A.
ZI Les Paluds
Avenue de Jouques – BP 1051
13781 Aubagne Cedex

Phone +33.442.845600
Fax +33.442.845619

Sartorius Stedim France SAS
ZI Les Paluds
Avenue de Jouques – CS 71058
13781 Aubagne Cedex

Phone +33.442.845600
Fax +33.442.846545

Austria

Sartorius Stedim Austria GmbH
Franzosengraben 12
A-1030 Vienna

Phone +43.1.7965763.18
Fax +43.1.796576344

Belgium

Sartorius Stedim Belgium N.V.
Leuvensesteenweg, 248/B
1800 Vilvoorde

Phone +32.2.756.06.80
Fax +32.2.756.06.81

Denmark

Sartorius Stedim Nordic A/S
Hoerskaetten 6D, 1.
DK-2630 Taastrup

Phone +45.7023.4400
Fax +45.4630.4030

Hungary

Sartorius Stedim Hungária Kft
Kagyló u. 5
2092 Budakeszi

Phone +36.23.457.227
Fax +36.23.457.147

Italy

Sartorius Stedim Italy S.p.A.
Via dell'Antella, 76/A
50012 Antella-Bagno a Ripoli (FI)

Phone +39.055.63.40.41
Fax +39.055.63.40.526

Netherlands

Sartorius Stedim Netherlands B.V.
Edisonbaan 24
3439 MN Nieuwegein

Phone +31.30.6025080
Fax +31.30.6025099

Poland

Sartorius Stedim Poland Sp. z o.o.
ul. Wrzesinska 70
62-025 Kostrzyn

Phone +48.61.647.38.40
Fax +48.61.879.25.04

Spain

Sartorius Stedim Spain SA
C/Isabel Colbrand 10,
Oficina 70
Polígono Industrial de Fuencarral
28050 Madrid

Phone +34.90.2110935
Fax +34.91.3589623

Switzerland

Sartorius Stedim Switzerland AG
Ringstr. 24 a
8317 Tagelswangen

Phone +41.52.354.36.36
Fax +41.52.354.36.46

U.K.

Sartorius Stedim UK Limited
Longmead Business Park
Blenheim Road, Epsom
Surrey KT19 9 QQ

Phone +44.1372.737159
Fax +44.1372.726171

America

USA

Sartorius Stedim North America Inc.
5 Orville Drive
Bohemia, NY 11716

Toll-Free +1.800.368.7178
Fax +1.631.254.4253

Sartorius Stedim SUS Inc.
1910 Mark Court
Concord, CA 94520

Phone +1.925.689.6650
Toll Free +1.800.914.6644
Fax +1.925.689.6988

Sartorius Stedim Systems Inc.
201 South Ingram Mill Road
Springfield, MO 65802

Phone +1.417.873.9636
Fax +1.417.873.9275

Argentina

Sartorius Argentina S.A.
Int. A. Avalos 4251
B1605ECS Munro
Buenos Aires

Phone +54.11.4721.0505
Fax +54.11.4762.2333

Brazil

Sartorius do Brasil Ltda
Av. Dom Pedro I, 241
Bairro Vila Pires
Santo André
São Paulo
Cep 09110-001

Phone +55.11.4451.6226
Fax +55.11.4451.4369

Mexico

Sartorius de México S.A. de C.V.
Circuito Circunvalación Poniente No. 149
Ciudad Satélite
53100 Naucalpan, Estado de México

Phone +52.5555.62.1102
Fax +52.5555.62.2942

Asia | Pacific

Australia

Sartorius Stedim Australia Pty. Ltd.
Unit 5, 7-11 Rodeo Drive
Dandenong South Vic 3175

Phone +61.3.8762.1800
Fax +61.3.8762.1828

China

Sartorius Stedim Beijing
Representative Office
No. 33, Yu'an Road,
Airport Industrial Zone B, Shunyi District
Beijing 101300

Phone +86.10.80426516
Fax +86.10.80426580

Sartorius Stedim Shanghai
Representative Office
Room 618, Tower 1, German Centre,
Shanghai, PRC., 201203

Phone +86.21.28986393
Fax +86.21.28986392.11

Sartorius Stedim Guangzhou Office
Room 704, Broadway Plaza,
No. 233-234 Dong Feng West Road
Guangzhou 510180

Phone +86.20.8351.7921
Fax +86.20.8351.7931

India

Sartorius Stedim India Pvt. Ltd.
#69/2-69/3, Jakkasandra
Kunigal Road, Nelamangala Tq
Bangalore – 562 123

Phone +91.80.4350.5361
Fax +91.80.4350.5253

Japan

Sartorius Stedim Japan K.K.
KY Building, 8-11
Kita Shinagawa 1-chome
Shinagawa-ku
Tokyo 140-0001

Phone +81.3.3740.5407
Fax +81.3.3740.5406

Malaysia

Sartorius Stedim Malaysia Sdn. Bhd.
Lot L3-E-3B, Enterprise 4
Technology Park Malaysia
Bukit Jalil
57000 Kuala Lumpur

Phone +60.3.8996.0622
Fax +60.3.8996.0755

Singapore

Sartorius Stedim Singapore Pte. Ltd.
1 Science Park Road,
The Capricorn, #05-08A,
Singapore Science Park 2
Singapore 117528

Phone +65.6872.3966
Fax +65.6778.2494